


GUIDED SURGERY

# KIT GSL GUIDED SINUS LIFT

Guided surgery simple and fast



Compact kit for simpler and more economical operation

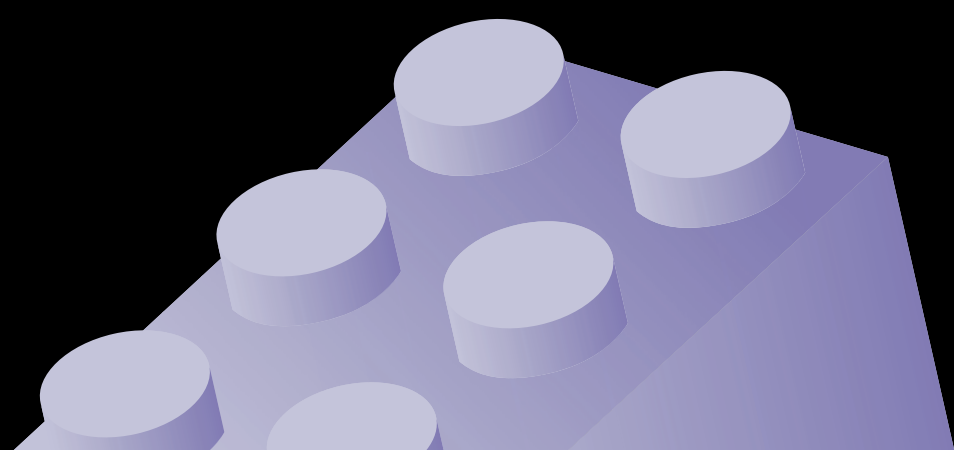

# PRODUCT DESCRIPTION

The **GSL Kit** is a system that allows the maxillary sinus to be augmented atraumatically and gradually to a height of 12 mm above its initial level, with a simple, quick and safe technique, performed in guided surgery and mixed technique. The bending of the cortical bone and the overcoming of the phase of elastic deformation until the bone breaks to allow the insertion of a reconstruction biomaterial and an implant are guaranteed by gradual and atraumatic series of steps of 1 mm each, using depth stops.

The **surgical instruments in GSL kit** are the following:

- 6 drills, with various profiles;
- 12 depth stops;
- 4 guided tips, different sizes
- 1 manual tool that has a straight working axis on one side and an "L" on the other that allows you to operate more easily even in cases of reduced opening;
- 1 instruments box with an autoclavable internal tray.

---

## Cylindrical Drills




## Chamfered Drills




## Rounded Drills



## Guided tips

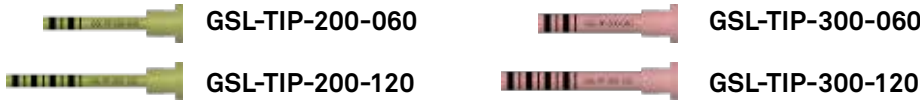

## Stops

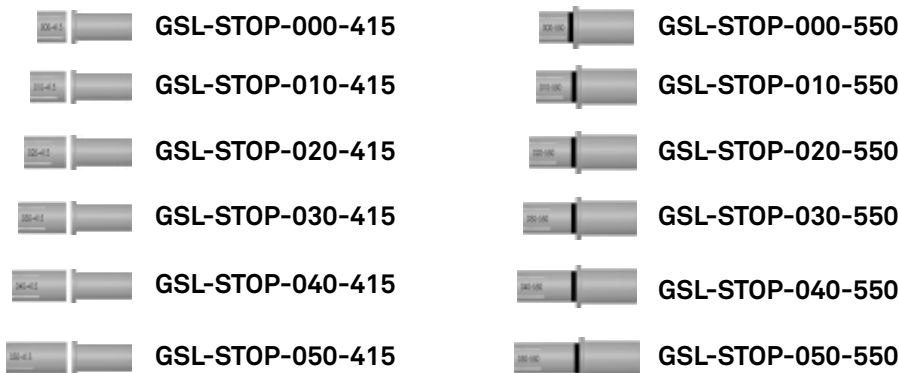

## Handle

